



AlcoHAWK[®] Slim

Owner's Manual

TABLE OF CONTENTS

INTRODUCTION	1
Overview	1
Basics of Breath Alcohol Testing	1
Effects of Alcohol	1
PREPARATION	3
Install the Batteries	3
Insert a New Mouthpiece	3
COMPONENTS DIAGRAM	4
OPERATION	5
READ THE RESULTS	6
ERROR MESSAGES	6
PRECAUTIONS	7
SPECIFICATIONS	8
VERIFICATION AND CALIBRATION	8
WARRANTY	9
DISCLAIMER	10

INTRODUCTION

Overview

The AlcoHAWK® Slim is a breath alcohol tester designed to test for the presence of alcohol in the breath. This product estimates blood alcohol concentration (BAC) in seconds. Please read this entire manual before using this product.

Basics of Breath Alcohol Testing

When a person consumes alcohol, the human body absorbs alcohol into the bloodstream. The alcohol in the bloodstream then exchanges with air deep in the lungs. The percent of alcohol in a person's blood (BAC) is proportional to the percent of alcohol in a person's breath (BrAC). Breath alcohol testers measure a person's deep-lung air alcohol level and convert it into BAC.

Effects of Alcohol

Variations among individuals and other factors, such as altitude and air temperature, affect the degree of intoxication. Some individuals may become intoxicated at low BAC levels. A low BAC on a breath alcohol tester does not mean that a person's reaction times can respond to any emergency encountered. Individuals should not drink alcohol and drive or use any machinery after consuming alcohol. There is no safe way to drive under the influence of alcohol. Even one drink can make you unsafe.

BAC	Dose-Specific Effects²
0.020-0.039% ³	Generally no loss of coordination. Slight euphoria or loss of shyness. Depressant effects are not apparent.
0.040-0.06%	General feeling of well-being, relaxation, lower inhibitions, and sensation of warmth. Euphoria. Some minor impairment of reasoning and memory. Lowering of caution. Driving skills may be impaired.
0.070-0.09%	Slight impairment of balance, speech, vision, reaction time, and hearing. Euphoria. Judgment and self-control are reduced. Caution, reason and memory are impaired. Driving skills are always impaired.
0.100-0.129%	Significant impairment of motor coordination and loss of good judgment. Speech may be slurred. Balance, vision, reaction time and hearing will be impaired. Euphoria.
0.130-0.159%	Gross motor impairment and lack of physical control. Blurred vision and major loss of balance. Euphoria reduced and dysphoria begins to appear.
0.160-0.249%	Dysphoria (anxiety, restlessness) predominates. Nausea may appear. The drinker has the appearance of a "sloppy drunk."
0.250%-0.299%	Needs assistance in walking. Total mental confusion. Dysphoria with nausea and some vomiting.
0.300%-0.399%	Loss of consciousness.
0.400% +	Onset of coma, possible death due to respiratory arrest. ⁴

² Bailey, William J., *Drug Use in American Society*, 3rd ed., Minneapolis: Burgess, 1993.

³ Some individuals may experience impairment at BAC levels at .03% or even lower.

⁴ Death can occur at a lower BAC level in some individuals.

PREPARATION

Install the Batteries

Apply a small amount of pressure on the indentation of the battery compartment cover located on the back of the product. Push down and away from the product to release the cover. Install 2 AA batteries and slide compartment cover back in place.




Opening and Inserting the Mouthpiece

Pull up on the folding mouthpiece located on the left side of the unit until the mouthpiece extends straight out from the unit. If testing multiple individuals, insert a new cover over the mouthpiece before each test.



COMPONENTS DIAGRAM

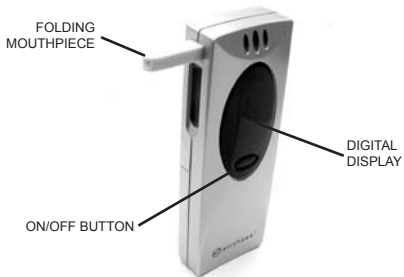

OPERATION



Wait 20 minutes after the test subject smokes, eats or drinks before testing. Failure to wait 20 minutes can give inaccurate BAC readings and damage the sensor.

- Step 1)** Place a fitted cover over the mouthpiece each time you test another individual.
- Step 2)** If it has been more than 24 hours since the device was last used, **prime the sensor** by blowing into the mouthpiece for 4 to 5 seconds **before** turning on the device.
- Step 3)** **PRESS** the black Power On/Off button located on the front of the unit. A single beep will sound and the device will count down from 60 as it prepares the sensor for testing.
- Note:** The device will restart a count down from 10 if the sensor is not ready within the first 60 second countdown.
- Step 4)** **WAIT** until the **RED LIGHT IS SOLID RED** and you hear a **SINGLE BEEP**.
- Note:** the device will not shut off until a test is performed OR until 15 seconds of non use elapses.
- Step 5)** Take a **DEEP BREATH AND BLOW STEADILY** (for at least 5 seconds) until a **BEEP** sounds to signal completion of the testing.

READ THE RESULTS

The displayed result is an estimate of the test subject's BAC. A breath alcohol tester cannot determine the exact BAC in the blood of the test subject. Correlation between BrAC and BAC depends on many factors, such as temperature and health conditions. Therefore, only use this device as a tool to screen for the presence of alcohol.

Many people consider a 0.02% or greater a positive test for alcohol in "zero tolerance" screening. Most states have 0.08% as the legal BAC limit for tasks such as operating a motor vehicle. Some states prohibit driving commercial vehicles or any vessel with a 0.04% BAC or higher. Many states prosecute individuals for driving under the influence at any level.

ERROR MESSAGES

Error Message	Cause and Course of Action
Sn or PL	There has been an error with the sensor. This can occur if the test subject does not wait at least 20 minutes after smoking, eating or drinking before using the unit or if the unit has not been used for an extended period of time. Press the Power On/Off button once to turn off the unit and repeat operation steps 1 – 5. If this message is displayed four or more consecutive times, the unit may need to be serviced.
rP	Repeat the test because the unit did not detect an even alcohol concentration in the breath sample. This can occur if the user varies the blowing force when performing a test or if the exhaust port is obstructed. Repeat operation steps 1 – 5 to perform another test.
H	The unit detected a very high alcohol concentration (above 0.40% BAC).
Lo	Low Battery. Replace AA batteries.

PRECAUTIONS

1. After drinking, smoking and eating, users should **wait at least 20 minutes** before testing.
2. Avoid testing in strong winds, in a closed room with a heavy amount of smoke, or where alcohol is present.
3. When the unit displays Lo, replace the AA batteries.
4. Do not blow cigarette smoke, food or liquid into the unit.
5. Only use the unit in a temperature range of 10 - 40°C or 50 - 104°F. Operation above or below this range may affect the results.
6. Avoid testing in the presence of any substances that contain methyl alcohol, isopropyl alcohol or acetone. These substances may interfere with the results of the test.
7. Conditions that increase the amount of ketones on the breath, such as diabetes and low caloric intake, may cause a false positive test.
8. BAC CHANGES OVER TIME DEPENDING ON MANY FACTORS. BAC AT ONE TIME CAN DIFFER FROM BAC A FEW MINUTES LATER.
9. DO NOT USE THIS PRODUCT WITH AN IGNITION INTERLOCK DEVICE OR WHILE INSIDE A VEHICLE.
10. YOU CANNOT USE THE RESULTS OF THIS PRODUCT IN COURT.
11. DO NOT USE THIS DEVICE TO DETERMINE IF YOU OR ANYONE ELSE CAN OPERATE ANY VEHICLE OR MACHINERY.
12. DO NOT LET YOURSELF OR ANYONE ELSE DRINK AND DRIVE. ALWAYS USE A DESIGNATED DRIVER WHEN YOU OR ANYONE ELSE DRINKS ALCOHOL.

SPECIFICATIONS

Dimensions	5" x 1.75" x .75" (12.7 cm x 4.45 cm x 1.91 cm)
Weight	1.8 oz (51 g) without batteries; 3.4 oz (96 g) with batteries
Housing	Shock resistant, molded plastic
Power Supply	Two (2) AA alkaline batteries
Battery Life	150 tests ¹
Sensor Technology	Sensitive Semiconductor Sensor
Blowing Time	5 seconds
Response Time	5 seconds ¹
Digital Display	2 digits (.00% BAC)
Operating Temperature	10 - 40°C (50 - 104°F)
Detection Range	.00% - .40% BAC
Air Sample	5 seconds
Operation	1 button
Warranty	1 Year Limited Warranty

¹ Estimates, results may vary.

VERIFICATION AND CALIBRATION

Calibration of the AlcoHAWK® Slim occurs during manufacturing using advanced alcohol simulation equipment. The manufacturer passing known alcohol concentrations through the sensor to set baseline values for testing. During testing, the unit compares users' breath samples to these baseline values. If the verification determines that the BAC is not consistent with actual BAC levels, a technician calibrates the unit. The accuracy of breath alcohol testers can fluctuate after six to twelve months, depending on frequency of use and operating conditions.

If the unit is providing inconsistent results, unusually high or low readings, or no readings at all, a professional needs to recalibrate the unit. **A technician should recalibrate this unit at least every six to twelve months. Monthly calibrations may be required if the unit is used daily, such as professional screening.**

For recalibration or support, please visit www.q3i.com or call 888-399-1687. During recalibration, a technician will perform a diagnostic check on the unit and reset the baseline sensor values.

ONE-YEAR LIMITED WARRANTY

Q3 Innovations, LLC (Q3I) warrants this unit to be free from defects in workmanship or material (excluding calibration) under normal use for ONE (1) YEAR from the date of retail purchase by the original purchaser ("Warranty Period"). **The warranty covers parts and labor on covered repairs; software updates if applicable; and ground freight back to the customer (USA only). Except as provided above or applicable law, Q3I provides this product "as is."**

If a defect arises and Q3I receives a valid claim within the Warranty Period, at its option, Q3I will either (1) repair the product, (2) exchange the product with a product that is new or which has been manufactured from new or serviceable used parts and is at least functionally equivalent to the original product, or (3) refund the purchase price of the product. When a refund is given, the return product to Q3I becomes Q3I's property. The warranty does not cover freight to Q3I. The warranty is void if anyone tampers with, maliciously damages, or physically abuses the unit.

This warranty is limited to the original consumer purchaser and is not transferable to, or enforceable by, any subsequent owner.

For warranty service during the Warranty Period, please follow the below instructions:

1. Call 888-399-1687 or email service@q3i.com to obtain an RMA number and shipping instructions. Return the unit postage paid, insured

and in suitable packaging.

2. For your protection, obtain a proof of delivery receipt. Shipping costs are your responsibility.
3. You must enclose the following with the unit:
 - a. The completed RMA form (available at www.q3i.com/support.php)
 - b. Your name, complete return address and written description of the problem. (No PO Box please.)
 - c. Your telephone number to reach you during normal business hours.
 - d. A copy of your dated sales receipt or invoice.

DISCLAIMER

Q3I makes no warranties, expressed or implied, as to the ability of this unit to determine whether, or to the extent to which, a test subject's mental or physical functionality, or judgment, is impaired by blood alcohol including whether the test subject is intoxicated, under any definition of that word. Q3I expressly disclaims any liability for incidental, special, or consequential damages of any nature under any legal theory. Any act or failure to act based upon a reading from this product is at the test subject's or user's own risk. Q3I assumes no responsibility for consequences to users or test subjects who use this product and are later under the influence of alcohol or to have had their judgment or any mental or bodily function impaired by alcohol. Correlation between BrAC and BAC depends on many variables, such as temperature and health conditions of the test subject. A safe or low BAC reading does not mean that the test subject's physical or mental performance or judgment (including reaction times) can respond to any emergency encountered. A breath alcohol tester cannot determine the exact concentration of alcohol in the blood of a test subject.

Q3I does not warrant that the operation of the product will be error-free. Q3I is not responsible for any consequences arising from failure to follow instructions or relating to the product's use. This warranty does not apply: (a) to consequences caused by accident, abuse, tampering, misuse, flood, fire, earthquake or other external causes; (b) to consequences caused by operating the product outside the permitted or intended uses described by Q3I; (c) to consequences caused by service performed by anyone who is not a representative of Q3I or a Q3I Authorized Service Provider; or (d) to a product or part that has been modified to alter functionality or capability without the written permission of Q3I.

TO THE EXTENT PERMITTED BY LAW, THIS WARRANTY AND THE REMEDIES SET FORTH ABOVE ARE EXCLUSIVE AND IN LIEU OF ALL OTHER WARRANTIES, REMEDIES AND CONDITIONS, WHETHER ORAL OR WRITTEN, STATUTORY, EXPRESS OR IMPLIED. AS PERMITTED BY APPLICABLE LAW, Q3I SPECIFICALLY DISCLAIMS ANY AND ALL STATUTORY OR IMPLIED WARRANTIES, INCLUDING, WITHOUT LIMITATION, WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE AND WARRANTIES AGAINST HIDDEN OR LATENT DEFECTS. IF Q3I CANNOT LAWFULLY DISCLAIM STATUTORY OR IMPLIED WARRANTIES, THEN TO THE EXTENT PERMITTED BY LAW ALL SUCH WARRANTIES ARE LIMITED TO THE DURATION OF THE EXPRESS WARRANTY AND TO THE REPAIR OR REPLACEMENT SERVICE AS DETERMINED BY Q3I IN ITS SOLE DISCRETION.

Q3I does not authorize any Q3I reseller, agent, or employee to make any modification, extension, or addition to this warranty. If anyone determines any terms to be illegal or unenforceable, the legality or enforceability of the remaining terms shall not be affected or impaired.

EXCEPT AS PROVIDED IN THIS WARRANTY AND TO THE MAXIMUM EXTENT PERMITTED BY LAW, Q3I IS NOT RESPONSIBLE FOR DIRECT, SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES UNDER ANY LEGAL THEORY, INCLUDING BUT NOT LIMITED TO LOSS OF USE; LOSS OF REVENUE OR INCOME; LOSS OF ACTUAL OR ANTICIPATED PROFITS (INCLUDING LOSS OF PROFITS FROM A CONTRACT); LOSS OF THE USE OF MONEY; LOSS OF ANTICI-


PATED SAVINGS; LOSS OF BUSINESS; LOSS OF OPPORTUNITY;
LOSS OF GOODWILL, OR LOSS OF REPUTATION.

Some countries, states and provinces do not allow the exclusion or limitation of incidental or consequential damages or allow limitations on how long an implied warranty or condition may last, so the above limitations or exclusions may not apply to you. This warranty gives you specific legal rights, and you may also have other rights that vary by country, state or province. This limited warranty is governed by and construed under the laws of the country where the product was purchased.

DO NOT DRINK ALCOHOL AND DRIVE

ALWAYS USE A DESIGNATED DRIVER
AFTER DRINKING ALCOHOL

BE SAFE AND RESPONSIBLE

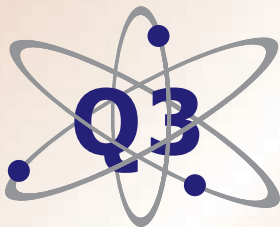

Q3 INNOVATIONS, LLC
Personal Safety & Monitoring Devices

2349 Jamestown Ave., Suite #4 • Independence, IA 50644

Tel: 888-399-1687 • Fax: 319-334-3421

Sales Email: sales@q3i.com • Service Email: service@q3i.com

www.q3i.com